

EddyCus® lab 2020A – Anisotropy Tester

P 2020A 23



Highlights

- ► Contact-free and real time
- ► Accurate single-point measurement
- Manual mapping guided by easy-tohandle software
- ► Measurement of encapsulated layers
- Characterization of multilayer materials upon request

Device Series

- ► Metal thickness (nm, µm)
- ► Sheet resistance (Ohm/sq)
- ► Emissivity
- ► Conductivity / resistivity (mOhm·cm)
- ► Electrical anisotropy (%)
- ▶ Weight (g/m²) and drying status (%)
- ▶ Permeability (H/m) Beta

SURAGUS GmbH Maria-Reiche-Strasse 1 01109 Dresden Germany

For further questions: +49 351 32 111 520

sales@suragus.com

Visit us at:
www.suragus.com
www.suragus.com/calculator
www.suragus.com/EddyCusLab2020

Applications

- ► Touch panel sensors (TPS)
- ► Printed electronics
- ▶ Wearable electronics
- ► Smart textiles
- ► Photovoltaics
- ► Smart / switchable films
- ► Medical surfaces and devices
- ► Biological sensors
- ► Aerospace, automotive, transportation
- ► Semiconductor and memory
- ► Energy storage

Materials

- ► Nanowire films
 - ► Conductive NW (Ag, Ni, Pt, Au)
 - ► Semiconductor NW (Si, SiC)
 - ► Magnetic NW (Fe₃O₄-AgNWs)
 - ► Multilayer NW (ZnO/AgNW/ZnO)
- ► Carbon Nano Tubes and Buds
- ► Fiber reinforced composites
- ► Metal meshes, smart meshes
- Anisotropic grain / domain materials
- Anisotropic effect / defect directions (cracks, line defects)


Engineered and Made in Germany






Anisotropy Term and Concept

- ► Electrical anisotropy refers to a difference in electrical resistance depending on the direction of current flow
- ▶ Wire and mesh structures can have anisotropic resistances
- ► Bulk materials with dominant directional characteristics / effects / defects can also have electrical anisotropy
- ► Anisotropy can be optimized to the layout of the contacts
- Anisotropy can save material and improve optical transparency to sheet resistance ratio

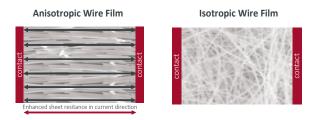

- ▶ Described by anisotropy direction and strength
- ▶ Both characteristics must be obtained at the same position
- ► The anisotropy strength is calculated using the lowest and highest resistance that align in perpendicular directions
- ► Inline deposition, e.g. slot die coating on moving web, tends to create lower resistances in machine direction "MD" and higher resistance in traversing direction "TD"
- ► Calculation as ratio of lowest and highest resistance

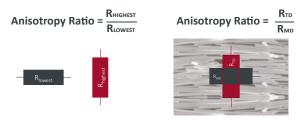

Device Characteristics

Measurement technology	Non-contact eddy current sensor with directed current induction
Substrates	Foils, glass, wafer, etc.
Substrate area	8 inch / 204 mm x 204 mm (open on three sides)
Max. sample thickness / sensor gap	3/5/10/25 mm (defined by the thickest sample)
Sheet resistance range	0.01 – 1,000 Ohm/sq; 1 to 5 % accuracy
Anisotropy range (TD/MD)	0.33 – 3 (larger upon request)
Device dimensions (w/h/d) / weight	11.4" x 5.5" x 17.5" / 290 mm x 140 mm x 445 mm / 10 kg
Further available features	Metal thickness, sheet resistance, emissivity, resistivity, weight and drying status and also permeability (beta) measurement

Device Control and Software

